

NU SKIN® RENUSPA iO™



Nu Skin RenuSpa iO delivers more of the benefits every body needs to feel their best. It's the integrated wellness and beauty device that gives you confidence as real as your results.

MEET NU SKIN RENUSPA iO

Say hello to a self-care win with the power of Nu Skin RenuSpa iO—our first fully integrated wellness and beauty system. The exclusive design features Adaptive Microcurrent Technology that comfortably adjusts to your skin in real time, plus an IoT-connected app, and perfectly paired Nu Skin® Conductive Gel keep you meeting and exceeding your goals.

Unleash the freedom to feel whole again with Nu Skin RenuSpa iO.

WHY YOU'LL LOVE IT

RESTORE BENEFITS

Through skin stimulation, helps:

- Depuff the skin's appearance.
- The skin feel energized.
- Provide a spa-like experience to lift mood, improve overall feelings of wellbeing, and help you look and feel better.

REVITALIZE BENEFITS

Through skin stimulation, helps:

- Visibly contour and reduces the appearance of dimples for a smoother, radiant look.
- Visibly tone and sculpt.
- Visibly firm and tighten for a more youthful appearance.

REFRESH BENEFITS

- Helps your body to feel refreshed and relaxed.
- Provides a smooth, comfortable massage while contributing to skin stimulation.

DEVICE FEATURES

- The surface topography and current provide a smooth, comfortable experience, while contributing to skin stimulation.

- The patent-pending Adaptive Microcurrent Technology continuously and comfortably adjusts to your skin in real time—for an optimized, just-for-you treatment, every time.
- The smooth, soft silicone sleeve makes cleaning simple and minimizes mess.

NU SKIN VERA® APP FEATURES

- A library of routines for all of your Restore, Revitalize, and Refresh needs.
- Guided Tutorials to help you every step of the way with result-driven instructions
- View live data insights through M-Power Insights which displays, in real-time, the power of Adaptive Microcurrent Technology.
- Delivers complementary product recommendations to help you reach your personal skincare goals.
- Designed to provide a personalized iO Score to help you stay motivated and reach your wellness and beauty goals.
- Developed to track your unique usage with fun, incentivizing badges.

CLINICAL RESULTS THAT PROVE IT

Forty-two subjects between 25–65 years of age were recruited to complete a Nu Skin RenuSpa iO regimen for 12 weeks.

The subjects used the RenuSpa iO device with Nu Skin Conductive Gel five times a week for five minutes per area (e.g., back of upper arm, back of upper thigh (buttocks), and abdomen). To assess visible skin depuffing, toning, and sculpting, thigh and arm circumferences were taken at baseline and weeks 2, 4, 8, and 12. Statistically significant reductions were noted throughout the 12 weeks of the study.

**Note: The use of Nu Skin RenuSpa iO does not contribute to weight loss.*


NU SKIN® RENUSPA iO™

Comparing use after 12 weeks to baseline, statistically relevant data shows the ability of the Nu Skin RenuSpa iO device, through skin stimulation, to visibly tone, hydrate, and reduce the appearance of dimpled skin on the upper thigh and abdomen.


Subjects also noted improvements in their skin's appearance within the following attributes for each area:

- Firmness
- Smoothness
- Contour


Week 12 Self-Perception—Skin Appearance

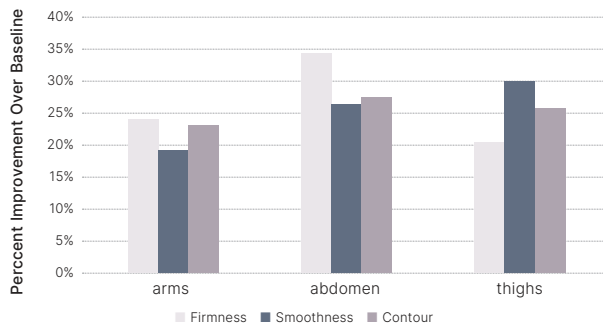


Expert Grader Hydration

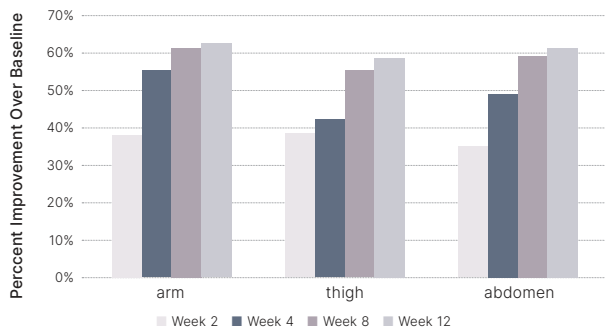

Expert Grader Tonicity

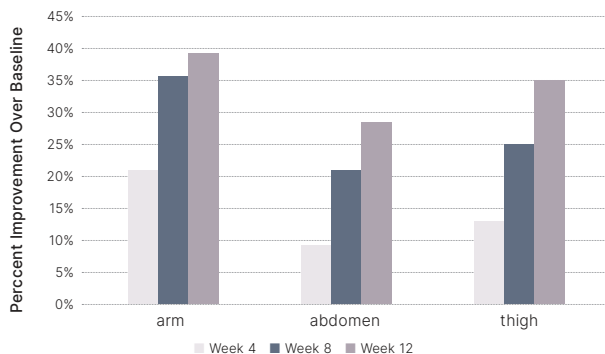

Expert Grader Appearance of Dimpled Skin

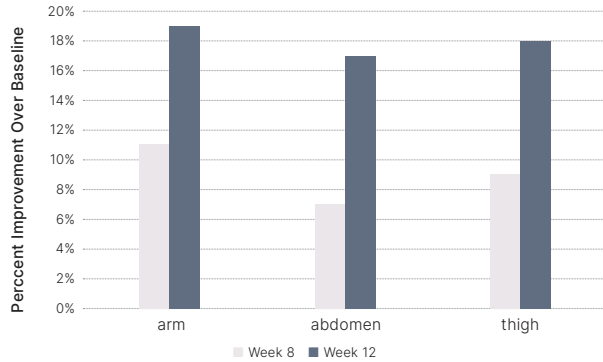

HOW TO USE IT

NU SKIN RENUSPA iO + NU SKIN CONDUCTIVE GEL

Use Nu Skin RenuSpa iO with Nu Skin® Conductive Gel on clean skin. Use Conductive Gel on clean skin. For best results, use the device 3 to 5 times a week per area. Do not use the device daily. Each session can last up to 5 minutes. See the User's Manual or connect to Nu Skin Vera for full usage instructions. Reapply gel as needed during the session to ensure smooth, continuous connectivity between the skin and device. Wipe off any remaining product with a damp cloth.

For best results, cleanse with ageLOC Dermatic Effects Polish as the first step, followed by ageLOC Dermatic Effects Body Contouring Lotion at the end of your routine.

WHAT POWERS IT

Conductive Gel is a refreshing formula used with your Nu Skin RenuSpa iO routines. It's formulated to work with the device to help you look and feel your best.

YOU MIGHT ALSO LIKE

- ageLOC® Dermatic Effects® Polish
- ageLOC® Dermatic Effects® Serum
- ageLOC® Dermatic Effects® Body Contouring Lotion
- Beauty Focus® MultiBeauty
- Beauty Focus® Collagen+
- ageLOC Meta®
- Epoch® IceDancer® Invigorating Leg Gel
- Sunright® Insta Glow

LEARN MORE ABOUT IT

How do you integrate wellness and beauty in one device?

More than ever, we're treating our bodies like an ecosystem and want holistic solutions that complement our personal health and evolving needs. Nu Skin RenuSpa iO does just that.

Through skin stimulation provided by the Adaptive Microcurrent

NU SKIN® RENUSPA iO™

Technology and its unique metallic nodes, RenuSpa iO provides wellness and beauty routines to supercharge your body inside and out. Plus, the spa-like experience helps to:

- Enhance wellbeing
- Energize the skin
- Uplift mood
- Refresh and relax the body
- Comfortably massage

Enhances skin appearances and visibly:

- Smooths
- Tones & sculpts
- Firms
- Increases skin radiance
- Depuffs skin

What are the key skin benefits I'll get from using Nu Skin RenuSpa iO?

The RenuSpa iO system is designed to deliver powerful wellness and beauty results through skin stimulation and its restoring, revitalizing, and refreshing benefits:

RESTORE—Boost your mood, depuff and energize your skin, and feel your best any time of day as Nu Skin RenuSpa iO gives you a spa day feeling of total-body wellbeing.

REVITALIZE—Help visibly firm, tighten, sculpt, and tone—while helping reduce the appearance of dimpled skin.

REFRESH—Comfortably massage your skin and help your body feel relaxed and refreshed.

What benefits will I enjoy because of Adaptive Microcurrent Technology?

The short answer is that it helps you maximize your wellness and beauty goals. But to be technical, the microcurrent continuously and comfortably adjusts to your skin 80+ times per second—for a personalized treatment every time.

What makes the Nu Skin RenuSpa iO system unique?

RenuSpa iO is a complete reimagination of the at-home spa experience. It utilizes a powerful Adaptive Microcurrent Technology that continuously adapts to the changes in your skin conductivity throughout each session to provide you with an optimal and comfortable experience. RenuSpa iO also incorporates one-of-a-kind, patent-pending metallic node design that helps provide a smooth, comfortable massage while contributing to skin stimulation. It's designed to perfectly fit the contours of your body and glide smoothly over your skin. Through skin stimulation, this system not only helps to contour and beautify the skin, but also helps depuff the skin's appearance and relax and refresh your body for an overall uplifting experience as you personalize your wellness and beauty journey.

What is the maximum amount of time recommended to use RenuSpa iO?

Use Nu Skin RenuSpa iO with Nu Skin Conductive Gel on clean skin. For best results, use the device 3 to 5 times a week for up to five minutes per area. Do not use the device daily.

Why are there different usage movements for RenuSpa iO?

Users looking for restoring benefits will focus on basic upward/inward movements of the entire area. When looking to revitalize your appearance, use movements that are directed up and/or inward and focus on more targeted areas, but also provide some gentle massaging of the area to help improve the appearance of the skin. During a Refresh routine, specific body areas are being targeted through gentle massage to help you feel overall relaxed and refreshed. Therefore, smaller, circular, upward/inward motions are desired.

Why do the metallic nodes have such a unique shape?

We designed the Nu Skin RenuSpa iO device's one-of-a-kind metallic nodes to help provide a smooth, comfortable massage while contributing to skin stimulation. Plus, we designed them to perfectly fit to the contours of your body and glide smoothly over your skin.

Do I have to connect with the app to use my Nu Skin RenuSpa iO?

If you don't connect, you'll still enjoy great wellness and beauty benefits. However, Nu Skin RenuSpa iO is created to connect, so you'll want to benefit from all the features developed exclusively for you, including personalized goal setting, guided tutorials, live data insights, and usage tracking.

How do I know which routine is right for me?

You know your wellness and beauty needs best, and there are plenty of routine options to choose from once you're in the Vera® app. Try them all to get all the amazing benefits we're sure your skin and body will love.

NU SKIN® RENUSPA iO™

How does Nu Skin RenuSpa iO support Nu Skin's commitment to sustainability?

As responsible stewards of our planet, we have proudly designed Nu Skin RenuSpa iO with a rechargeable battery and to be more water resistant, increasing the longevity of the device. Please check your local device recycling capabilities and recycle the carton where possible. By purchasing Nu Skin RenuSpa iO, you're joining us in taking an active stride to ensure the future is bright for both the planet and each of us.

Please visit www.nuskin.com/sustainability for more information and details on how to recycle.

WHAT'S IN IT

NU SKIN CONDUCTIVE GEL

Water (Aqua), Glycerin, Pentylene Glycol, Carbomer, Sodium Hydroxide, Chlorphenesin.



Scan the QR code to connect with Vera®

For more information, please consult the Nu Skin® RenuSpa iO User Manual

